

**Bring Big Cleaning  
Power to Small Spaces**

The Hamo® LS-2000  
Washer/Disinfector



## Documented Cleaning in a Small Space

STERIS's longstanding reputation for unparalleled product performance and service, combined with Hamo's tradition of innovative product design, has now resulted in a stronger, more comprehensive range of washer/disinfectors that can meet all of your cleaning needs. This provides:

- > State-of-the-art systems that balance durable, precision engineering with effective chemical cleaning.
- > Units that are manufactured in compliance with the latest U.S. and European standards.
- > Total solutions that deliver higher value and performance to lower overall costs, and environmental impact.
- > Meets prEN/ISO 15883-1/2 and HTM 2030 standards.

## The Hamo LS-2000 Washer/Disinfecter

Efficient, powerful and highly configurable, the Hamo LS-2000 Washer/Disinfecter brings big equipment cleaning capability to small spaces. The LS-2000 offers the largest processing capacity for a machine of its size.

### The Flexibility to Match Your Unique Space and Tasks


Small but highly productive, the LS-2000 provides you with flexibility to get the most of out of your limited space.

- > Unique racks and washing manifolds accommodate nearly every type of reusable medical device.
- > Set at the appropriate height for transfer carts that make loading and unloading simple.
- > A wide range of validated STERIS detergents ensure effective and thorough reprocessing of delicate instruments.
- > A glass window enables observation of load processing.

### Superior Throughput and Cleaning Efficiency

The LS-2000 offers best-in-class capacity for its size, so you can process more critical instruments and equipment faster.

- > The capacity-expanding, 5-level washing manifold can process up to 10 DIN trays per cycle.
- > Rapid, dirty-to-clean pass-through functionality eliminates time consuming pre-cleaning steps without compromising productivity.
- > Low water and detergent consumption yields more efficient cycles compared to larger units with a smaller capacity.
- > The high-powered drying system speeds drying time and helps reduce environmental impact.

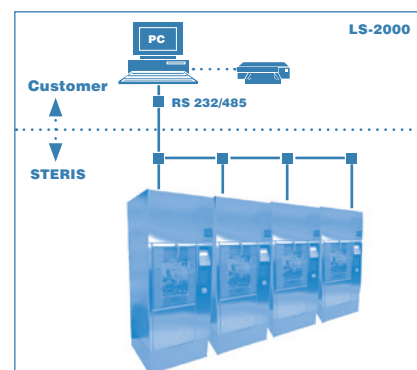

*IDMS (Intelligent Data Management System) or PDVS (Process Data Visualisation and Validation Management) options are available.*





### Bring Big Cleaning Power to Small Spaces

The Hamo LS-2000 Washer/Disinfector, with its big cleaning capacity and space-saving size, brings many benefits to your institution.

It's quiet operation and ample processing flexibility makes the LS-2000 a perfect choice for your cleaning requirements.

- > Operating rooms
- > GI suites
- > Outpatient clinics
- > Ambulatory surgery centres
- > Reprocessing

### High-Performance Cleaning, Predictable Results

Every cycle run through the Hamo LS-2000 Washer/Disinfector meets the water and time parameters required for thermal disinfection, so you can count on consistent outcomes cycle after cycle.

- > A highly efficient HEPA filter ensures effective drying without airborne contamination.
- > A professional analysis of your facility's water can optimise detergent input, enabling consistent cycles and outcomes.
- > Up to 99 programmable cycles offer unprecedented control over your cleaning program.

### Built-in Data Management Provides Increased Compliance, Connectivity, and Control

Robust data management and control software extends the benefits of the LS-2000 well beyond its small footprint. The powerful Intelligent Data Microprocessor Control software enables managers to pre-set standard cleaning processes to eliminate user-error and help ensure compliance with regulatory standards.


- > All systems can link via RS-232 (via SCADA protocol) to any hospital-wide tracking system improving traceability.
- > The optional Process Data Validation System (PDVS) software validates each cycle to help ensure compliance with prEN-15883 standards.
- > The optional IDMS visualisation software allows you to remotely monitor air and water temperatures—providing the peace of mind that the process is on-track.
- > The Hamo LS-2000 Endoclean option enables safe, efficacious processing of flexible endoscopes. (Available in selected markets)



## Worldwide Service and Support

From systems planning to industry-leading technical service, all STERIS washer/disinfectors are backed by an uncompromising commitment to our customers.

- > STERIS trained professionals are dedicated to ensuring that every STERIS system delivers maximum performance.
- > 24-hour engineering and technical service support is just a telephone call away. (Available in selected markets)
- > A full range of service contracts can help extend the life of your equipment.
- > Experienced system planners can consult with your facility to design the most efficient CSSD process.

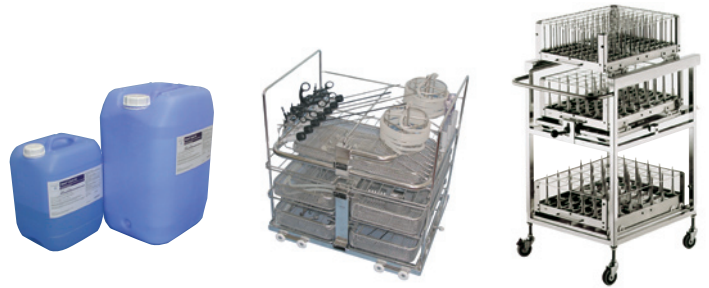

*For detergents, racks and conveyors, please ask for specific brochures.*

## Effective Solutions for Hospital Decontamination

	Maximum Rack Dimension			Control Systems		Applications		
	Width mm (inch)	Height mm (inch)	Depth mm (inch)	Micro-Processor	PC Siemens	Endoclean Option	Multi-Chamber	Ultrasonic
<b>Ecoline® LM-15</b>	420 (16.5)	285 (11.2)	460 (18.1)	✓				
<b>Ecoline® LM-25</b>	550 (21.7)	570 (22.4)	500 (19.7)	✓				
<b>Ecoline® BM-10</b>	500 (19.7)	360 (14.2)	550 (21.7)	✓				
<b>Hamo® LS-950</b>	480 (18.9)	500 (19.7)	480 (18.9)	✓		✓		
<b>Hamo® LS-2000</b>	500 (19.7)	630 (24.8)	620 (24.4)	✓		✓		
<b>Hamo® T-21</b>	640 (25.2)	640 (25.2)	640 (25.2)	✓		✓		✓
<b>Hamo® T-21L</b>	640 (25.2)	640 (25.2)	850 (33.5)	✓				✓
<b>Hamo® T-840</b>	640 (25.2)	640 (25.2)	640 (25.2)		✓		✓	✓
<b>Hamo® T-840L</b>	640 (25.2)	640 (25.2)	850 (31.5)		✓		✓	✓
<b>Hamo® DW Cart Washer</b>	1000 (39.4)	2000 (78.8)	1600 (63) 2700 (106.3)	✓				

STERIS is committed to helping healthcare professionals implement highly effective, environmentally friendly washer/disinfectant systems that optimise productivity, reliability, and safety.

To learn more about the Hamo LS-2000 Washer/Disinfectant, or any of our technologies to prevent infection and contamination, visit [www.steris.com](http://www.steris.com), or call your local in-country office.

STERIS recommends the following validated cleaning chemistries for use with the Hamo LS-2000 Washer/Disinfectant. (Available in selected markets)

- > Hamo® Liquid 52 Neutral Enzymatic Detergent
- > Hamo® Liquid 36 Biocompatible Drying Aid
- > Hamo® Liquid 50 Alkaline Detergent
- > Hamo® Liquid 55 Alkaline Detergent
- > Hamo® Liquid 84 Acidic Neutraliser

Made in Switzerland by  
STERIS AG, Bielstrasse 76  
CH-2542 Pieterlen/Biel  
[www.hamo.com](http://www.hamo.com)

Publication ID #M2526GB.2005-07, Rev. B  
GPSI Printed 07/2005, 2,500  
©2005 by STERIS Corporation  
All rights reserved. Printed in USA

Technologies to Prevent Infection and Contamination™

### STERIS OFFICES WORLDWIDE

Canada	800 661 3937
France	33 1 55 176000
Germany	49 2233 6999 0
Italy	39 02 2130341
Japan	81 78 321 2271
Korea	82 2 517 1517
Latin America	305 442 8202
Singapore	65 68 41 7677
Spain	34 91 658 5920
Sweden	46 152 228 30
Switzerland	41 32 376 02 00

# STERIS®



[www.steris.com](http://www.steris.com)

### Global Headquarters

STERIS Corporation  
5960 Heisley Road  
Mentor, OH 44060-1834 ■ USA  
440-354-2600 ■ 800-548-4873

### European Headquarters

STERIS Limited  
STERIS House  
Jays Close  
Viables  
Basingstoke  
Hampshire RG22 4AX  
+44 1256 840400  
UK Freephone 0800 252609